

# Hematopoietic Progenitor Cell (CD34<sup>+</sup>) Total RNA

## human

Cell type-specific total RNA for gene cloning and gene expression analysis      Order no. 130-093-167

### Contents

1. Description
  - 1.1 Background information
  - 1.2 Applications
  - 1.3 Reagent and instrument requirements
2. Product quality
  - 2.1 Blood donation
  - 2.2 Cell preparation
  - 2.3 Cell type-specific total RNA preparation and purity
3. Examples
4. References

### 1. Description

<b>Components</b>	0.2 µg Hematopoietic Progenitor Cell (CD34 <sup>+</sup> ) Total RNA, human (lyophilized).
<b>Capacity</b>	0.2 µg total RNA.
<b>Source</b>	Healthy human blood donors.
<b>Storage</b>	Dissolve cell type-specific total RNA in 10 µL sterile, RNase-free distilled water and prepare aliquots. Store dissolved total RNA and lyophilized total RNA at -70 °C.

#### 1.1 Background information

CD34 is expressed on hematopoietic stem cells and all late progenitor cells such as pre-B cells, colony forming unit erythrocytes (CFU-E) and colony forming unit granulocytes (CFU-G). It is not expressed on pluripotent stem cells (CD133<sup>+</sup>CD34<sup>-</sup>).

Hematopoietic Progenitor Cell (CD34<sup>+</sup>) Total RNA has been isolated from highly purified human CD34<sup>+</sup> hematopoietic progenitor cells.

#### 1.2 Applications

The RNA is suitable for gene expression profiling, cDNA library generation, RT-PCR analysis, cloning, and characterization of progenitor cell-specific genes.

#### 1.3 Reagent and instrument requirements

- Sterile, RNase-free distilled water.

### 2. Product quality

#### 2.1 Blood donation

Donors were consenting, healthy, normal adults aged between 18 and 68 years. The donation procedure was performed following the guidelines of the German Medical Association and the regulations of German and EU law. Donors were negative for HIV, hepatitis B and C, and syphilis. Donors undergoing a course of medication were excluded.

#### 2.2 Cell preparation

Peripheral blood mononuclear cells (PBMCs) were prepared from buffy coats, collected from multiple donors, by Ficoll™ density gradient centrifugation. CD34<sup>+</sup> progenitor cells were isolated from the PBMCs by positive magnetic selection using the CD34 MicroBead Kit, human (# 130-046-702). Purities always exceeded 90% CD34<sup>+</sup> cells.

#### 2.3 Cell type-specific total RNA preparation and purity



Pooled, purified cells from multiple donors were lysed with RNA lysis buffer and total RNA was extracted by using silica-membrane technology. The RNA was treated with RNase-free DNase I to remove residual contamination with genomic DNA. After treatment DNase was removed and inactivated. RNA purity was determined by capillary electrophoresis with an Agilent Bioanalyzer and consistently showed RNA Integrity Number (RIN) values<sup>1,2</sup> of over 8.5.

Table 1: Donor, cell, and RNA quality parameters

Parameter	Value
Donor-age distribution	18–68 years
Donor gender	50% male, 50% female
Viral serology of donors	Negative for HIV, HBV HCV, syphilis
Donor medication	Negative
Cell purity	> 90% CD34 <sup>+</sup>
RNA Integrity Number (RIN) <sup>1,2</sup>	> 8.5

### 3. Examples

The Hypoxanthine-guanine phosphoribosyltransferase (HPRT) gene is a housekeeping gene that is normally expressed at a low frequency of 1–10 mRNA copies per cell,<sup>3</sup> therefore, it provides a suitable reference for RNA quality. To determine HPRT mRNA content 1 ng, 10 ng, and 100 ng of Hematopoietic Progenitor Cell (CD34<sup>+</sup>) Total RNA were analyzed by real-time RT-PCR using intron-spanning primers.

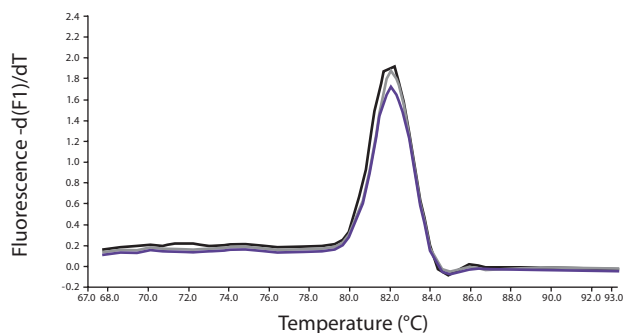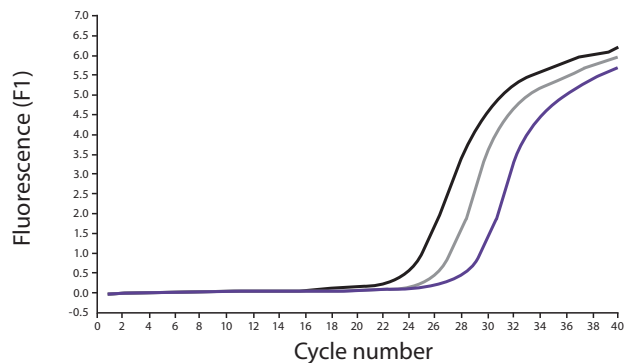

### 4. References

1. Schroeder, A. *et al.* (2006) The RIN: an RNA integrity number for assigning integrity values to RNA measurements. *BMC Molecular Biology* 2006, 7: 3
2. Imbeaud, S. *et al.* (2005) Towards standardization of RNA quality assessment using user-independent classifiers of microcapillary electrophoresis traces. *Nucl. Acids Res.* 2005 33: e56.
3. Steen, A. M. *et al.* (1990) Levels of hypoxanthine phosphoribosyltransferase RNA in human cells. *Exp. Cell Res.* 186: 236–244.

All protocols and data sheets are available at [www.miltenyibiotec.com](http://www.miltenyibiotec.com).

#### Warranty

The products sold hereunder are warranted only to be free from defects in workmanship and material at the time of delivery to the customer. Miltenyi Biotec GmbH makes no warranty or representation, either expressed or implied, with respect to the fitness of a product for a particular purpose. There are no warranties, expressed or implied, which extend beyond the technical specifications of the products. Miltenyi Biotec GmbH's liability is limited to either replacement of the products or refund of the purchase price. Miltenyi Biotec GmbH is not liable for any property damage, personal injury or economic loss caused by the product.

MACS is a registered trademark of Miltenyi Biotec GmbH.

Ficoll is a trademark of GE Healthcare companies.

Copyright © 2009 Miltenyi Biotec GmbH. All rights reserved.