

## Contents

1. Description
  - 1.1 Background information
  - 1.2 Applications
  - 1.3 Recommended antibody dilution
  - 1.4 Reagent requirements
2. General protocol for immunofluorescent staining
3. Examples of immunofluorescent staining with Mouse IgG2a isotype control antibodies

## 1. Description

**Components** 1 mL monoclonal Mouse IgG2a isotype control antibodies, conjugated to various dyes.

|            |             |
|------------|-------------|
| FITC       | 130-091-837 |
| PE         | 130-091-835 |
| APC        | 130-091-836 |
| VioBlue®   | 130-094-671 |
| VioGreen™  | 130-096-923 |
| PerCP      | 130-094-967 |
| PE-Vio770™ | 130-096-638 |
| APC-Vio770 | 130-096-637 |

The antibody concentrations are 22 µg/mL (FITC-conjugate), 22 µg/mL (PE-conjugate), 11 µg/mL (APC-conjugate), 55 µg/mL (VioBlue-conjugate), 55 µg/mL (VioGreen-conjugate), 55 µg/mL (PerCP-conjugate), 11 µg/mL (PE-Vio770-conjugate), or 22 µg/mL (APC-Vio770-conjugate).

**Clone** S43.10 (isotype: mouse IgG2a).

**Capacity** 100 tests or up to 10<sup>9</sup> total cells.

**Product format** Antibodies are supplied in buffer containing stabilizer and 0.05% sodium azide.

**Storage** Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

### 1.1 Background information

The Mouse IgG2a isotype control antibody clone S43.10 is specific for the hapten NP (4-hydroxy-3-nitro-phenyl) acetyl. This hapten is not expressed on cells or cell lines. Therefore, the antibody clone S43.10 can be used as a negative control, to distinguish specific from non-specific binding of mouse IgG2a fluorochrome-conjugated antibodies to human cells, for example, via Fc receptors, or due to interactions of the fluorochrome with the cell surface.

### 1.2 Applications

- Mouse IgG2a isotype control antibodies are suitable for assessing the level of background staining for flow cytometric cell analysis of cells.

### 1.3 Recommended antibody dilution

The recommended antibody dilution for all Mouse IgG2a conjugates is **1:11 for up to 10<sup>7</sup> cells/100 µL** of buffer for labeling of cells and analysis by flow cytometry. The antibody is suited for staining of formaldehyd-fixed cells.

Cells should be stained prior to fixation, if formaldehyde is used as a fixative.

### 1.4 Reagent requirements

- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS® BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Keep buffer cold (2–8 °C).

▲ **Note:** EDTA can be replaced by other supplements such as anticoagulant citrate dextrose formula-A (ACD-A) or citrate phosphate dextrose (CPD). BSA can be replaced by other proteins such as human serum albumin, human serum, or fetal bovine serum (FBS). Buffers or media containing Ca<sup>2+</sup> or Mg<sup>2+</sup> are not recommended for use.

- (Optional) Propidium Iodide Solution (# 130-093-233) or 7-AAD for flow cytometric exclusion of dead cells without fixation.
- (Optional) Fixation and Dead Cell Discrimination Kit (# 130-091-163) for cell fixation and flow cytometric exclusion of dead cells.









## 2. General protocol for immunofluorescent staining

▲ Volumes given below are for up to 10<sup>7</sup> nucleated cells. When working with fewer than 10<sup>7</sup> cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 2×10<sup>7</sup> nucleated cells, use twice the volume of all indicated reagent volumes and total volumes).

1. Determine cell number.
2. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
3. Resuspend up to 10<sup>7</sup> nucleated cells per 100 µL of buffer.
4. Add 10 µL of the Mouse IgG2a isotype control antibody.
5. Mix well and incubate for 10 minutes in the dark in the refrigerator (2–8 °C).
  - ▲ **Note:** Higher temperatures and/or longer incubation times may lead to non-specific cell labeling. Working on ice requires increased incubation times.
6. Wash cells by adding 1–2 mL of buffer and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
7. Resuspend cell pellet in a suitable amount of buffer for analysis by flow cytometry or fluorescence microscopy.

### 3. Examples of immunofluorescent staining with Mouse IgG2a isotype control antibodies

Human peripheral blood lymphocytes were stained with Mouse IgG2a isotype control antibodies conjugated to FITC (A), PE (B), APC (C), VioBlue (D), VioGreen (E), PerCP (F), PE-Vio770 (G), or APC-Vio770 (H) and analyzed by flow cytometry using the MACSQuant® Analyzer. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence. Cells were gated on viable lymphocytes.

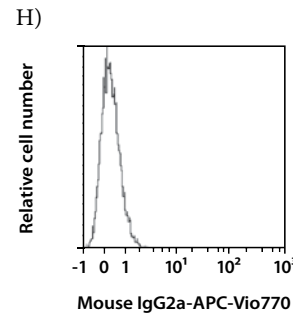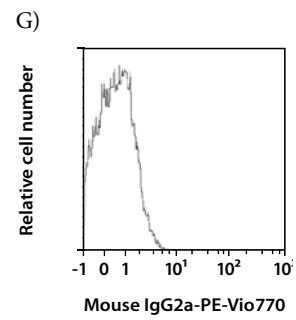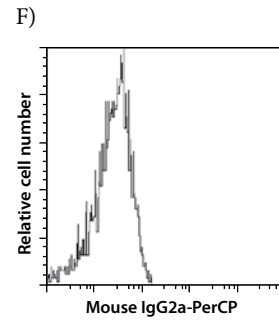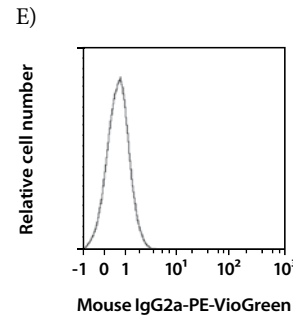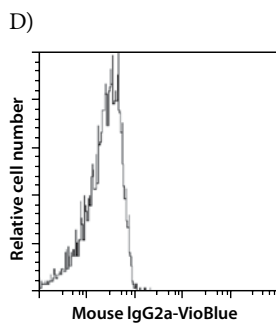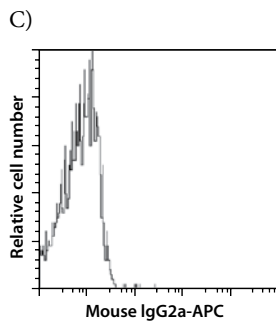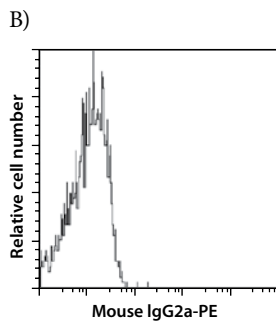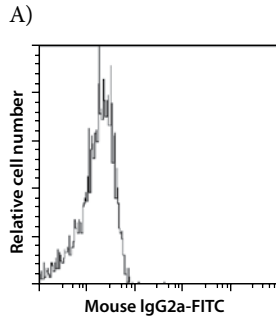

All protocols and data sheets are available at [www.miltenyibiotec.com](http://www.miltenyibiotec.com).

### **Warnings**

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

### **Warranty**

The products sold hereunder are warranted only to be free from defects in workmanship and material at the time of delivery to the customer. Miltenyi Biotec GmbH makes no warranty or representation, either expressed or implied, with respect to the fitness of a product for a particular purpose. There are no warranties, expressed or implied, which extend beyond the technical specifications of the products. Miltenyi Biotec GmbH's liability is limited to either replacement of the products or refund of the purchase price. Miltenyi Biotec GmbH is not liable for any property damage, personal injury or economic loss caused by the product.

autoMACS, MACS, MACSQuant, Vio, and VioBlue are registered trademarks and VioGreen and Vio770 are trademarks of Miltenyi Biotec GmbH.

Copyright © 2011 Miltenyi Biotec GmbH. All rights reserved.