



Miltenyi Biotec

Anti-ACSA-2 MicroBead Kit

mouse

Order no.
Order no.

130-097-678
130-097-679

Contents

1. Description
 - 1.1 Principle of the MACS® Separation
 - 1.2 Background information
 - 1.3 Applications
 - 1.4 Reagent and instrument requirements
2. Protocol
 - 2.1 Sample preparation
 - 2.2 Magnetic labeling
 - 2.3 Magnetic separation
3. Example of a separation using the Anti-ACSA-2 MicroBead Kit
4. References

Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

1. Description

This product is for research use only.

Components	50 µL Anti-ACSA-2 MicroBeads, mouse (# 130-097-679) or 1 mL Anti-ACSA-2 MicroBeads, mouse (# 130-097-678); MicroBeads conjugated to monoclonal anti-ACSA-2 antibodies (isotype: rat IgG2b). 50 µL FcR Blocking Reagent, mouse (# 130-097-679) or 1 mL FcR Blocking Reagent, mouse (# 130-097-678)
Capacity	For 5×10 ⁷ total cells, up to 5 separations (# 130-097-679) or for 10 ⁹ total cells, up to 100 separations (# 130-097-678).
Product format	All reagents are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

1.1 Principle of the MACS® Separation

First, Fc receptors are blocked with FcR Blocking Reagent, mouse. Then, the ACSA-2⁺ cells are magnetically labeled with Anti-ACSA-2 MicroBeads. The cell suspension is loaded onto a MACS® Column, which is placed in the magnetic field of a MACS Separator. The magnetically labeled ACSA-2⁺ cells are retained within the column. The unlabeled cells run through; this cell fraction is thus depleted of ACSA-2⁺ cells. After removing the column from the magnetic field, the magnetically retained ACSA-2⁺ cells can be eluted as the positively selected cell fraction. To increase the purity, the positively selected fraction containing the ACSA-2⁺ cells can be separated over a second column.

1.2 Background information

Anti-ACSA-2 MicroBeads (ACSA-2: astrocyte cell surface antigen-2) have been developed for the isolation of astrocytes from cell suspensions of mouse neural tissue based on the expression of the ACSA-2 antigen.

The ACSA-2 antigen is expressed specifically on astrocytes in a pattern similar to GLAST. Therefore, ACSA-2 is a specific marker for astrocytes in the central nervous system (CNS). The percentage of ACSA-2⁺ astrocytes differ according to the mouse age and the brain region used for cell isolation.

Anti-ACSA-2 MicroBeads allow for the isolation of astrocytes from all mouse stages.

1.3 Applications

- Positive selection or depletion of cells expressing the mouse ACSA-2 antigen.
- Isolation or depletion of ACSA-2⁺ astrocytes from dissociated mouse brain tissue.

1.4 Reagent and instrument requirements

- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, and 0.5% bovine serum albumin (BSA) by diluting MACS® BSA Stock Solution (# 130-091-376) 1:20 with PBS. Keep buffer cold (2–8 °C). Degas buffer before use, as air bubbles could block the column.
 - ▲ Note: BSA can be replaced by other proteins such as mouse serum albumin, mouse serum, or fetal bovine serum (FBS).
- The Neural Tissue Dissociation Kit (P) (# 130-092-628) is recommended for the generation of single-cell suspensions of neural cells from mouse brain tissue from mice up to P7.
 - ▲ Note: If additional cell surface epitopes of interest are papain sensitive, the Neural Tissue Dissociation Kit (T) (# 130-093-231) can be used.
- Adult Brain Dissociation Kit, mouse and rat (# 130-107-677) is recommended for the generation of single-cell suspensions of neural cells from mouse brain tissue derived from mice older than P7.

- Pre-Separation Filters (70 μm) (# 130-095-823) to remove cell clumps.
- MACS Columns and MACS Separators: ACSA-2⁺ cells can be enriched by using MS or LS Columns or depleted with the use of LD Columns. Positive selection or depletion can also be performed by using the autoMACS[®] Pro Separator.

Column	Max. number of labeled cells	Max. number of total cells	Separator
Positive selection			
MS	10 ⁷	2×10 ⁷	MiniMACS, OctoMACS, VarioMACS, SuperMACS II
LS	2×10 ⁷	4×10 ⁷	MidiMACS, QuadroMACS, VarioMACS, SuperMACS II
Depletion			
LD	1.5×10 ⁷	3×10 ⁷	MidiMACS, QuadroMACS, VarioMACS, SuperMACS II
Positive selection or depletion			
autoMACS	5×10 ⁷	10 ⁸	autoMACS Pro

▲ **Note:** Column adapters are required to insert certain columns into the VarioMACS™ or SuperMACS™ II Separators. For details refer to the respective MACS Separator data sheet.

- (Optional) Fluorochrome-conjugated Anti-ACSA-2 antibodies for flow cytometric analysis, e.g., Anti-ACSA-2-PE (# 130-097-674). For more information about antibodies refer to www.miltenyibiotec.com/antibodies.
- (Optional) Propidium Iodide Solution (# 130-093-233) or 7-AAD for flow cytometric exclusion of dead cells.
- (Optional) MACS Neuro Medium (# 130-093-570) and MACS NeuroBrew[®]-21 (# 130-093-566) for cultivation.
- (Optional) gentleMACS™ Dissociator (# 130-093-235), gentleMACS Octo Dissociator (# 130-095-937), gentleMACS Octo Dissociator with Heaters (# 130-096-427), and gentleMACS C Tubes (# 130-093-237, # 130-096-334).

2. Protocol

2.1 Sample preparation

For the preparation of single-cell suspensions from neural tissues refer to the data sheet of the Neural Tissue Dissociation Kit (P) (# 130-092-628) or Adult Brain Dissociation Kit (# 130-107-677), which can be used in combination with the gentleMACS Dissociator (# 130-092-235).

▲ **Note:** If additional cell surface epitopes of interest are papain sensitive, the Neural Tissue Dissociation Kit (T) (# 130-093-231) can be used.

A special protocol for the isolation of astrocytes “Isolation of astrocytes from adult mouse brain” is available at our website www.miltenyibiotec.com/130-097-678 in the library tab.



2.2 Magnetic labeling

▲ Work fast, keep cells cold, and use pre-cooled solutions. This will prevent capping of antibodies on the cell surface and non-specific cell labeling.

▲ Volumes for magnetic labeling given below are for up to 10⁷ total cells. For best performance it is recommended to use at least 5×10⁶ cells. When working with fewer than 10⁷ cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 2×10⁷ total cells, use twice the volume of all indicated reagent volumes and total volumes).

▲ For optimal performance it is important to obtain a single-cell suspension before magnetic labeling. Pass cells through 70 μm nylon mesh (Pre-Separation Filters (70 μm), # 130-095-823) to remove cell clumps which may clog the column. Moisten filter with buffer before use.

▲ The recommended incubation temperature is 2–8 °C. Higher temperatures and/or longer incubation times may lead to non-specific cell labeling. Working on ice may require increased incubation times.

1. Determine cell number.
2. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
3. Resuspend cell pellet in 80 μL of buffer per 10⁷ total cells.
4. Add 10 μL of FcR Blocking Reagent per 10⁷ total cells.
5. Mix well. Do not vortex. Incubate for 10 minutes in the refrigerator (2–8 °C).
6. Add 10 μL of Anti-ACSA-2 MicroBeads per 10⁷ total cells.
7. Mix well. Do not vortex. Incubate for 15 minutes in the refrigerator (2–8 °C).
8. Wash cells by adding 1–2 mL of buffer per 10⁷ cells and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
9. Resuspend up to 10⁷ cells in 500 μL of buffer.
10. Proceed to magnetic separation (2.3 or 2.4).

▲ **Note:** For higher cell numbers, scale up buffer volume accordingly.



2.3 Magnetic separation

▲ Choose an appropriate MACS® Column and MACS Separator according to the number of total cells and the number of ACSA-2⁺ cells. For details refer to the table in section 1.4.

▲ Always wait until the column reservoir is empty before proceeding to the next step.

Magnetic separation with MS or LS Columns

1. Place column in the magnetic field of a suitable MACS Separator. For details refer to the respective MACS Column data sheet.

2. Prepare column by rinsing with the appropriate amount of buffer:

MS: 500 µL LS: 3 mL

3. Apply cell suspension onto the column. Collect flow-through containing unlabeled cells.

4. Wash column with the appropriate amount of buffer. Collect unlabeled cells that pass through and combine with the flow-through from step 3.

MS: 3×500 µL LS: 3×3 mL

▲ **Note:** Perform washing steps by adding buffer aliquots as soon as the column reservoir is empty.

5. Remove column from the separator and place it on a suitable collection tube.

6. Pipette the appropriate amount of buffer onto the column. Immediately flush out the magnetically labeled cells by firmly pushing the plunger into the column.

MS: 1 mL LS: 5 mL

7. (Optional) To increase the purity of ACSA-2⁺ cells, it is recommended to enrich the positive fraction over a second MS or LS Column. Repeat the magnetic separation procedure as described in steps 1 to 6 by using a new column.

▲ **Note:** Elution of the cells from the column after the separation should be performed with cell culture medium if cells are to be taken directly into culture, otherwise elute with buffer as before.

▲ **Note:** Keep handling time of cells in PBS/BSA buffer to a minimum.

Depletion with LD Columns

1. Place LD Column in the magnetic field of a suitable MACS Separator. For details refer to the LD Column data sheet.
2. Prepare column by rinsing with 2 mL of buffer.
3. Apply cell suspension onto the column.
4. Collect unlabeled cells that pass through and wash column with 2×1 mL of buffer. Collect total flow-through; this is the unlabeled cell fraction. Perform washing steps by adding buffer two times. Only add new buffer when the column reservoir is empty.

Magnetic separation with the autoMACS® Pro Separator

▲ Refer to the respective user manual for instructions on how to use the autoMACS® Pro Separator.

▲ Buffers used for operating the autoMACS Pro Separator should have a temperature of ≥10 °C.

1. Prepare and prime the instrument.
2. Apply tube containing the sample and provide tubes for collecting the labeled and unlabeled cell fractions. Place sample tube in row A of the tube rack and the fraction collection tubes in rows B and C.
3. For a standard separation choose one of the following programs:

Positive selection: Possel


Collect positive fraction in row C of the tube rack.

Depletion: Depletes

Collect negative fraction in row B of the tube rack.

3. Example of a separation using the Anti-ACSA-2 MicroBead Kit

ACSA-2⁺ cells were isolated from P4 CD-1 mouse brain tissue using the Neural Tissue Dissociation Kit (P), the gentleMACS Dissociator, FcR Blocking Reagent, mouse, Anti-ACSA-2 MicroBeads, a MiniMACS™ Separator, and an MS Column. Cells were fluorescently stained with Anti-ACSA-2-PE and analyzed by flow cytometry using the MACSQuant® Analyzer. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.

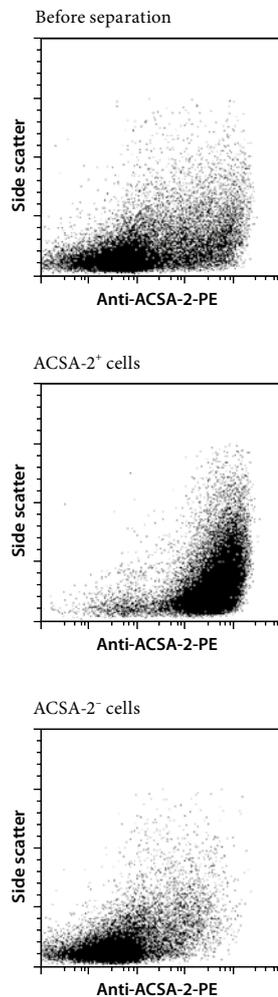

Refer to www.miltenyibiotec.com for all data sheets and protocols.

4. References

1. Feldmann, M. *et al.* (2014) Isolating astrocytes and neurons sequentially from postnatal murine brains with a magnetic cell separation technique. *Journal of Biological Methods* 1(2): e11.
2. Sharma, K. *et al.* (2015) Cell type- and brain region-resolved mouse brain proteome. *Nature neuroscience* 18(12): 1819–1831.
3. Holt, L. M. & Olsen, M. L. (2016) Novel applications of magnetic cell sorting to analyze cell-type specific gene and protein expression in the central nervous system. *PLoS one* 11(2): e0150290.
4. Kantzer, C. G. *et al.* (2017) Anti-ACSA-2 defines a novel monoclonal antibody for prospective isolation of living neonatal and adult astrocytes. *Glia* 65(6): 990–1004.
5. Batiuk, M. Y. *et al.* (2017) An immunoaffinity-based method for isolating ultrapure adult astrocytes based on ATP1B2 targeting by the ACSA-2 antibody. *J. Biol. Chem.* 292(21): 8874–8891.

Warranty

The products sold hereunder are warranted only to be free from defects in workmanship and material at the time of delivery to the customer. Miltenyi Biotec GmbH makes no warranty or representation, either expressed or implied, with respect to the fitness of a product for a particular purpose. There are no warranties, expressed or implied, which extend beyond the technical specifications of the products. Miltenyi Biotec GmbH's liability is limited to either replacement of the products or refund of the purchase price. Miltenyi Biotec GmbH is not liable for any property damage, personal injury or economic loss caused by the product.

autoMACS, gentleMACS, MACS, MACSmix, MACSQuant, MidiMACS, MiniMACS, MultiMACS, NeuroBrew, OctoMACS, QuadroMACS, SuperMACS, and VarioMACS are either registered trademarks or trademarks of Miltenyi Biotec GmbH and/or its affiliates in various countries worldwide.

CD-1 is a registered trademark of Charles River Laboratories International, Inc.

Copyright © 2017 Miltenyi Biotec GmbH and/or its affiliates. All rights reserved.